

**Bexen cardio**  
Ready for life



Reanibex  
**100**

**Feel  
Cardioprotected**

2 options:

- Semiautomatic
- Fully- Automatic

Adult and pediatric patient use


Alerts about the battery and pads status

CPR with metronome and feedback

Configured according to ERC/AHA guidelines

# Reanibex 100

Ready for life

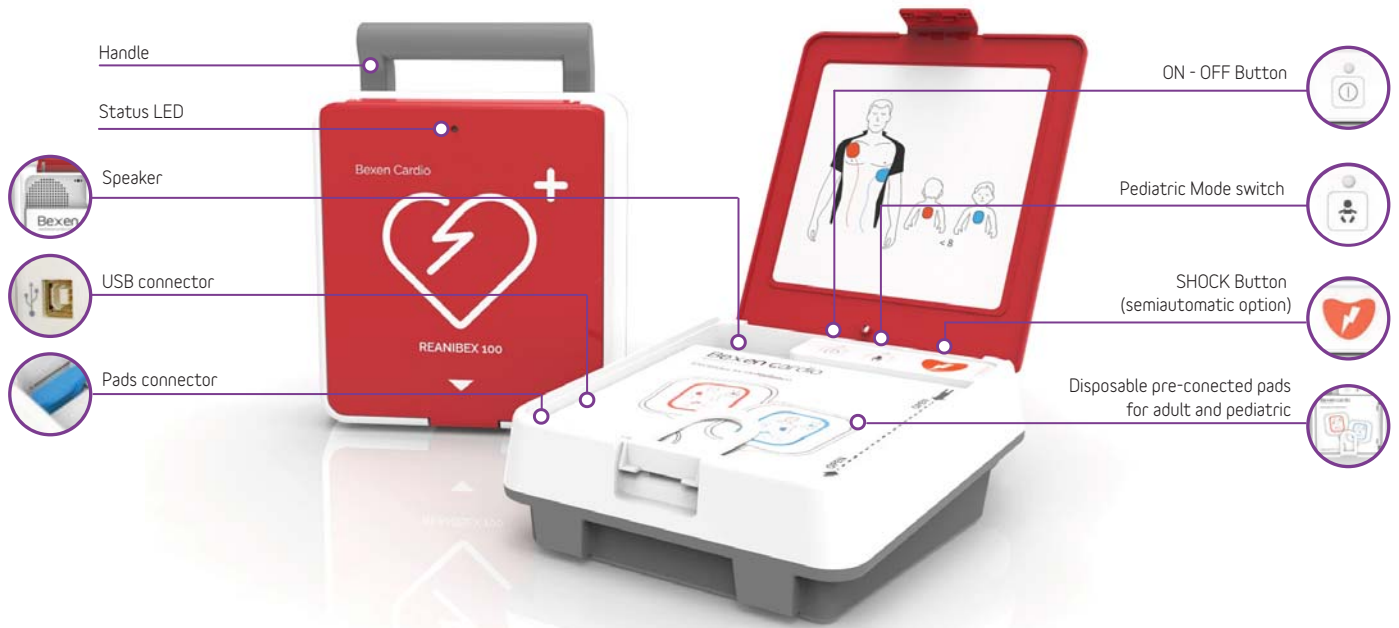

## FEATURES



### CLEAR GUIDANCE

Audible prompts to guide users through the complete operation sequence. Volume adapted to background noise.



### READY TO USE

The status LED flashes to indicate that the REANIBEX 100 is ready to be used. If a problem is detected the status LED is off.



### SEMI- OR FULLY-AUTOMATIC

In the Semi-automatic version the rescuer has to push the shock button to deliver the shock. In the fully automatic version the shock is delivered automatically when needed.



### FASTER SHOCK

Automatically switched on when the cover is opened, pre-connected pads and precharge for a faster shock delivery.



### SELF CHECK

Automatic daily selftest to ensure device and accessories continuous readiness.



### CPR SUPPORT

Metronome and rate feedback in real time for effective chest compressions.

## APPLICATIONS



### REANIBEX DATA NET

Web application for monitoring AED devices. Alerts generation about the status of the devices and their accessories, equipment configuration, software updates and the possibility to display test results and episodes (episodes optional).



### REANIBEX NFC READER

NFC application for Android devices to check the status of the REANIBEX 100 and its accessories without PC connection.



### REANIBEX USB SYNC

PC application to synchronize the AED with the Reanibex Data Net via USB to configure the device and to download the information stored in the REANIBEX 100.

**Bexen cardio**

Osatu S. Coop.  
Edificio Zearrekobuelta  
Subida Areitio, 5  
48260 Ermua (Spain)

Tel: +34 943 170 220  
Fax: + 34 943 170 227  
info@bexencardio.com  
www.bexencardio.com



Rev: 01/2021